

Instructions for Use

Product's name: Straps



Savaria-Dent Kft
H-9700 Szombathely
Dr. Szabolcs Z. u. 5/A.
HUNGARY
Tel.: +36 94/505-840



Applicationv

Function: A class I medical device that when connected to an orthodontic appliance, serves as a support for pulling.

Max. 2.5kg safety module release force load capacity per side.

Indications

To be used during orthodontic treatments, based on the doctor's recommendation.

Contraindications


Allergic reactions while wearing the headcap during treatments.

Material: 100% cotton

Washing under 40°C with general detergent used for cotton.

The product cannot be used again once there is damage in the material, or along the seam.

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

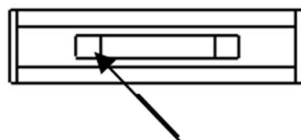

The safety module's connection

Version: 3.0 – 01.07.2021

Straps Product list

Ref. code	Trade name	Ref.code	Trade name
SD 48160	Neckpad + strap, wide ear	SD 48167	Wangenpolster, black
SD 48161	Strap jeans, narrow ear	SD 48170	Neckpad + strap, narrow ear
SD 48162	Strap, wide ear	SD 48171	Neckpad + strap, wide ear
SD 48162-L	Strap, wide ear L	SD 48173	Neckpad B + strap P, narrow ear
SD 48162-XL	Strap, wide ear XL	SD 48174	Neckpad B + strap Y, narrow ear
SD 48163	Strap, narrow ear	SD 48181	Jeans bunt, wide ear
SD 48163-L	Strap, narrow ear L	SD 48182	Jeans bunt, narrow ear
SD 48163-XL	Strap, narrow ear XL	SD 48186	Jeans bunt, jeans, narrow ear
SD 48164	Neckpad, buttonhole	SD 48187	Jeans bunt, jeans, wide ear
SD 48166	Wangenpolster, blue		