

## A Concentrated Experience

GC Orthodontics has been a distributor of self-ligating brackets from Japan-based TOMY Inc. since 2012.

The latest addition to the product line is EXPERIENCE™ C mini, offering all the advantages of its predecessor, EXPERIENCE™ C, but in a remarkably smaller size providing an incomparable aesthetic solution to patients.

EXPERIENCE™ C mini is manufactured by injection (Ceramic Injection Molded) from innovative ceramic materials (micro-particles). This new technology enables the development of smaller brackets, with a high resistance to impact and enhanced patient comfort.

EXPERIENCE™ C mini features a CNC milled slot for perfect torque expression and improved archwire engagement, along with a Rhodium-treated interactive clip for enhanced aesthetics. The chamfered slot facilitates archwire insertion, helping to eliminate archwire binding and notching. In addition, EXPERIENCE™ C mini's anatomical mechanical base follows the contours of the tooth to ensure a secure bond.



## Self Ligation Features

- Increased patient comfort
- Reduced chair time
- Fewer office visits necessary
- Shortened treatment time
- Improved patient hygiene
- Pleasing aesthetics

**SL** Product

**Hi** Aesthetic

# EXPERIENCE™

## mini ceramic

## Mini Self-Ligating Ceramic Bracket



Passive



Interactive



Active



**Experts in Orthodontics**

GC Orthodontics Europe GmbH


[www.gcorthodontics.eu](http://www.gcorthodontics.eu)

# EXPERIENCE mini ceramic

**Smooth swept tie-wings**  
engineered to minimize occlusal interference

**Multiple clip opening methods**  
for easy clip opening without clip deformation

**Ceramic injection molding bracket**  
for uniformity, strength & precision

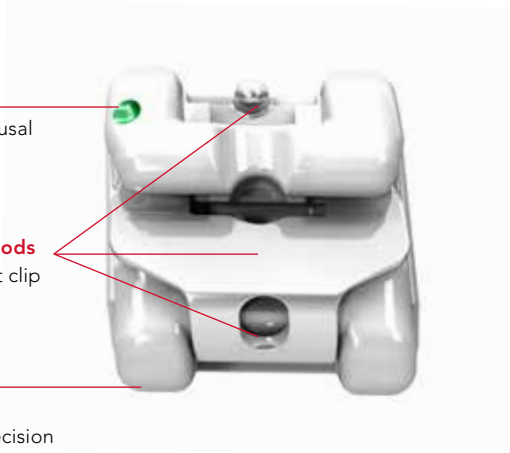

**CNC Milled slot**

**Triple chamfered slot**  
for easy wire engagement and to eliminate arch wire binding & notching

**Torque in Base**



**Compound contoured pad**  
ensures an anatomically correct fit for each tooth

**Mechanical lock**  
for maximum bond strength & the ability to reposition, if necessary. ID marking

## Bracket & Wire Combination

Pair EXPERIENCE with our Initalloy, Bio-Active, and Bio-Edge thermal NiTi arch wires, for increased patient comfort and treatment efficiencies. Rhodium arch wires also available for a completely discrete orthodontic bracket system.



The EXPERIENCE collection is manufactured by TOMY Inc., a leading manufacturer of self-ligating brackets since the year 2000. EXPERIENCE is available in both a standard and mini size; metal, rhodium, and ceramic material; and a variety of prescription options.

## EXPERIENCE mini ceramic Advantages

This ceramic self-ligating twin bracket delivers outstanding aesthetics and performance. The bracket features a rhodium coated clip that blends into the ceramic body of the bracket creating a beautiful appearance. The base is a 100% mechanical lock ensuring both reliable bond strength and easy debonding.

99-12ET-1300 · Rev00-06/2019



Since 1921  
Towards Century of Health

**GC Orthodontics Europe GmbH**  
Harkortstraße 2 · 58339 Breckerfeld · Germany  
Tel +49.2338.801.888 · info.gco.germany@gc.dental