




Features

- 15 to 20 % smaller
- Better patients comfort
- Céramique transparente, pure polycristalline alumina
- Thin Particles
- Perfect tooth mimic
- No bracket coloration
- Easy debonding notch

Hi Aesthetic

Chíc. mini

Mini ceramic bracket
Aesthetic Revolution

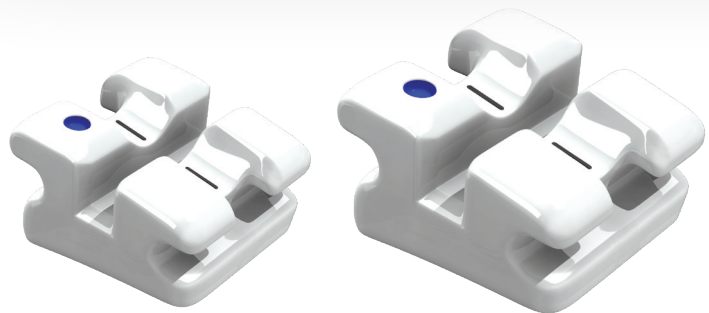

GC


Experts in Orthodontics

GC Orthodontics Europe GmbH

www.gorthodontics.eu




Chic. mini

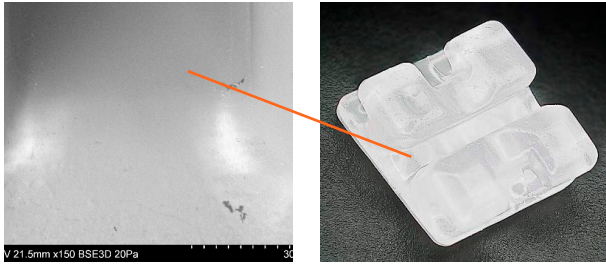


Low friction

The chamfered parts of the slot are injected in one piece, there is no setback due to the mold joint, the risk of buttressing will be less significant in severe crowding. Additionally, the archwire slide smoothly thanks to a unique slot's polishing process.

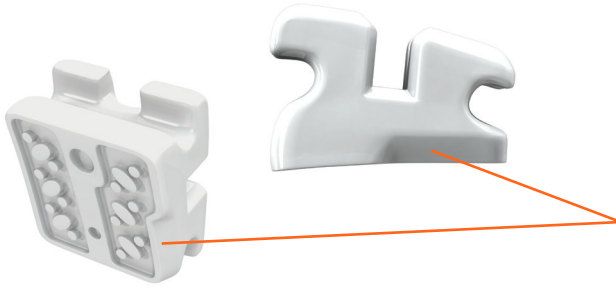

Easy Debonding

With its preformed notch on the base, Chic. mini is debonded safely without altering the bond strength as effectively as on a normal size bracket.

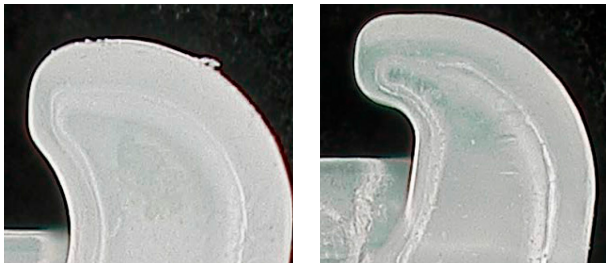

New Designed Hook

This new hook shape ensure a perfect chain or elastic retention.



Mini Size

This new design, on the lower incisors, increases the inter-bracket space, thus facilitating the insertion of the archwire at the first appointment, in cases of severe crowding.

Brackets combination

You can use Chic. mini on the upper with Legend Rhodium brackets on the lower. Combined with Rhodium archwires for a complete aesthetic result.



Chic. mini est fabriqué par TOMY Inc. in Yumoto, Japan. Orthodontists around the world have used and adored TOMY brackets since their first bracket system launched in December of 1961. As of today, TOMY has grown to be one of the most respected bracket manufacturers in the world.



99-15ET-1300 · Rev00 - 07.2022



GC Orthodontics Europe GmbH
Harkortstraße 2 · 58339 Breckerfeld · Germany
Tel +49.2338.801.888 · info.gco.germany@gc.dental